

A black and white photograph of a man in a wheelchair, wearing a tank top and shorts, holding a basketball high above his head with his right hand. The background is dark, and the lighting highlights the man's physique and the wheelchair's wheels. A bright yellow curved shape is at the bottom right.

QUICKIE
Live without limits.

ALL COURT

ALL COURT

Designed with the help of Sunrise Medical's elite wheelchair basketball athletes, the Quickie All Court® wheelchair establishes a new standard in performance and style. Using our patented center-of-mass adjustment feature, the All Court adjusts to match each user's unique center-of-gravity and seat-angle requirements - all without swapping a single part.

This, along with a new series of options, makes the Quickie All Court® the perfect chair for both the recreational and the elite athlete.



Available in Aluminium and Titanium versions

ALL COURT
Aluminium



The new casters ensure that chair speed and maneuverability are not simply absorbed in the wheels.

Extreme camber options of 16° and 20° offer superior stability and effortless turning.






Anti-tip Brace : The 5th wheel anti-tip brace is a new option featured on the Quickie All Court®. This adjustable brace, located between the backrest towel bar and the integrated anti-tip, eliminates the rearward flex in the seat frame caused when the player leans backwards. With this option, the player no longer has to adjust to the flex in the seat when leaning back to shoot over a defender.

The Quickie All Court® features single (standard) or dual (optional) integrated fifth-wheel anti-tip(s). These anti-tip caster(s) allow players to move their weight directly over the rear wheels, providing a tremendous advantage in maneuverability and power.

Our patented center-of-mass adjustment system provides a maximum horizontal (center-of-gravity) adjustment-range of 3.5" and up to 5" of vertical seat-squeeze (or seat dump) adjustment. This combination of adjustment-systems allows users to align their weight over the push-wheels and optimize their seat angles without making changes in fork size, caster size, caster-stem size, or fork-hole location.

ALL COURT
Titanium

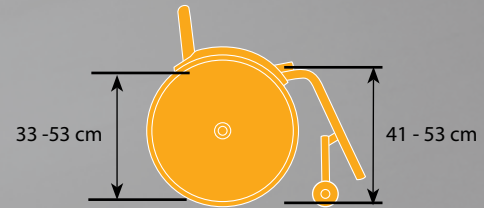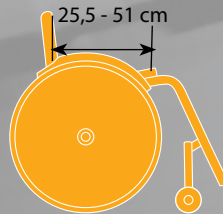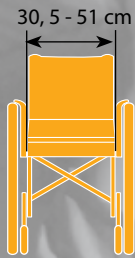

ALL COURT

ALL COURT TI

	ALL COURT	ALL COURT TI
Weight	from 10,88 kg	from 9,75 kg
Maximum user weight	113 kg	113 kg
Caster options	3" Micro	3" Micro
Rear wheel options	Spoke wheels and Spinergy sports wheels: 24" - 26"	Spoke wheels and Spinergy sports wheels: 24" - 26"
Axle options	Stainless steel, optional titanium or quad-release axles	Stainless steel, optional titanium or quad-release axles
Camber options	16° or 20°	16° or 20°

For further information on the full specification, options and accessories, please refer to the order form.

All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.com or alternatively is available on request in large text