

JAY

J3 BACK

- ≡ Now more than ever, we've got your back.



JAY J3 BACK

≡ Never before have the creative minds at JAY® spent more time studying the subtleties of fitting the back to the chair and the user to the back. The result: The J3 Back.

The most comprehensive combination of adjustable hardware, individual shell sizes and positioning capabilities ever assembled.

Three innovative hardware options: JayMount, LITE Mount and 4 Point hardware


Engineered to mount on crowded backposts ranging from 1 1/8" to 3/4", JAYMount hardware virtually eliminates incompatibility with multiple mounting locations. All three dimensions of height, angle and depth can be adjusted while the user remains in the chair. In addition to these great features, the LITE Mount hardware is aluminum and extremely lightweight. For increased stability, the 4 Point hardware can support more force to the back.

A multitude of shapes and sizes:

"We've got your back". Providing the correct size back for each individual is critical to preventing problems down the road. For this reason, the J3 back offers virtually every combination of depth, height and width, ranging from 12" to 26", that an adult user may require. This also includes the Heavy Duty shell, which combined with the 4 Point hardware has a weight capacity of 225 kg.



- Multiple angle adjustments without losing seat depth

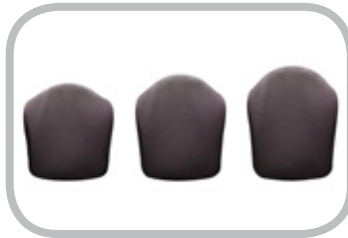

Spine-Align Posture Conforming System:

The Spine Align System is simple to use and creates the individual shape necessary to best support the user.

Result: Optimal postural support, pressure distribution and comfort.

Angle Adjustable Swing Away Laterals:

Accommodate a wider variety of positioning needs.



SPECIFICATIONS



JAY J3 Back	
Weight	0.77 kg and up
Product Width	12" – 26"
Product Depth	2.2" – 6"
Product Height	6.5" – 24"
Weight Capacity	136 kg. or 226 kg over 22" wide
Cover	Microclimatic 3DX™ Spacer Fabric